


Suction System Cleaning for your Sinus Chemical Ratio's

Daily –


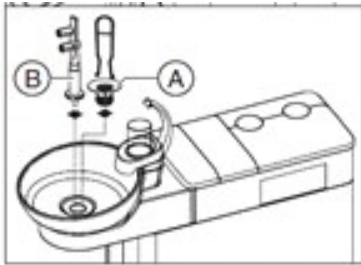


Alpro-Jet D or Orotol Plus - cleaning should be carried out at the end of each day *

Instructions



1. Remove the solids trap (A)
2. Empty the solids in to a suitable waste container
-  **Do not tip the solid in to the spittoon bowl with the solids trap removed**
3. Insert the cleaning adapter (B) as far as it will go

(*once a week Alpro-Jet D or Orotol Plus is replaced with Alpro-Jet W or MD 555. When Weekly cleaning is done Daily cleaning is NOT required)

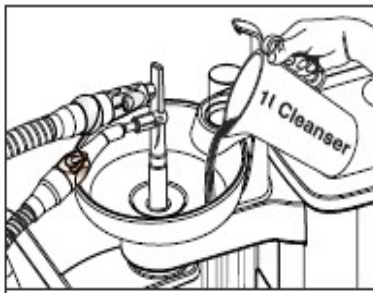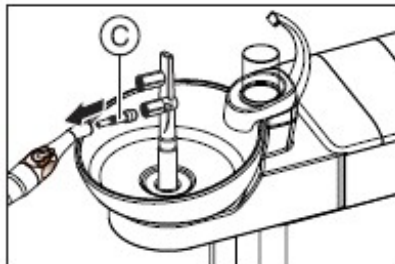

4. Prepare 1 litre of Alpro-Jet D or Orotol Plus cleaning solution in a separate container according to the manufactures instructions and mix it thoroughly.
- 50 mls of Alpro-Jet D : 950 mls of water
or
- 20 mls of Orotol Plus : 980 mls of water
5. Attached the HVE & LVE hoses to the cleaning adaptor (C)
6. Open the suction valves
7. Pour the cleaning mixture in to the spittoon bowl, the solution will be sucked through the hoses for approximately 4 mins. After 4 mins turn off the suction valves, lift the cleaning adaptor out and allow any excess solution to be drained away.
8. Return the cleaning adaptor back in to the spittoon bowl with the suction hoses attached and repeat the process but with 1 litre of fresh water.
9. After 4 mins turn replace the suction hoses back in to their holders and allow any excess water to be flushed down the spittoon.
10. Wipe out the spittoon bowl.

Weekly -

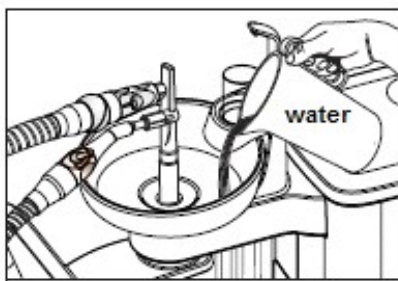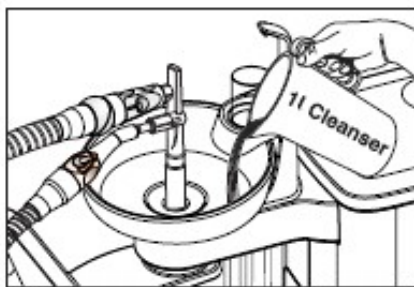

Alpro-Jet W or MD 555 - One day a week - at lunchtime

1. Follow the Alpro-Jet D / Orotol Plus instructions 1 – 3.
2. Prepare 1 litre of Alpro-Jet W or MD 555 cleaning solution in a separate container according to the manufactures instructions and mix it thoroughly.
- 50 mls of Alpro-Jet W : 950 mls of water
or
- 100 mls of MD 555 : 900 mls of water
3. Follow Alpro-Jet D / Orotol Plus Instructions 5 & 6.
4. Pour the Alpro-Jet W or MD 555 mixture in to the spittoon bowl and allow 2/3rd of the solution to be sucked up.
5. Disconnect the suction hoses and replace them back in their hangers
6. Allow the solution to soak in the hoses for 30 minutes, the remaining solution can be flushed down the spittoon bowl
7. After 30 mins reconnect the HVE & LVE hoses back on to the cleaning adaptor, flush the lines with approximately 2 litres of fresh water.
8. Wipe out the spittoon bowl

(Daily cleaning is NOT required on the day weekly cleaning is carried out)