

Vertex®

Articulating System

Works the way you work.



The Vertex Tray.

Saves time.

- Lets you create removable dies without traditional pin indexing.
- Ribs in the tray wall make it easy to guide the dies back into place and prevent mesial-distal movement.
- The tray spine, which becomes an integral part of the poured stone base, is engineered to prevent vertical and lateral displacement of the die, much like a pin.
- Preparation of the die and model can be completed in a single pour with grinding and drilling.



Provides Accuracy.

- The precisely engineered ribbed walls guide and stabilize the die while the spine provides a firm, positive seat.
- Each Vertex Tray is identical to every other Vertex Tray. You can simply transfer the die and model to a new tray and re-articulate, should it ever become necessary to adjust the bite.



Offers the Perfect Combination.

- Use the new Vertex Tray with the renowned Vertex Articulator to address the vast majority of articulations in your laboratory with speed, accuracy and economy.

The Vertex Articulator.

Offers Simplicity With Accuracy.

- The unique, patented design allows a full range of lateral and protrusive movements while maintaining accurate centric occlusion.
- Easy separation of the upper and lower models provides easy access to your work.

Delivers Even More Versatility.

- You can use it with your current pinning system or combine it with the Vertex Tray.
- The Vertex Articulator will streamline your production process, reduce the amount of material you use and enhance the overall quality of your work.

Makes A Presentation That Says "Professional."



- The Vertex Articulator and the Vertex Tray will give your cases the clean professional look that your dentists will appreciate.
- You can customize the mold so that every case carries your lab name.



Vertex Articulating System

Introductory Kit.....PC 100098

This kit contains 10 articulators (slot type), 10 articulators (flex slot type), 2 rubber base molds (posterior and medium-full), 6 trays and spines, 1 adhesive, cyanoacrylate (20g), 1 adhesive accelerator (2 oz.), 1 release block and 1 plastic stop rod.

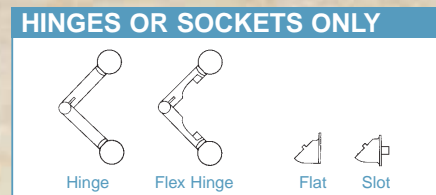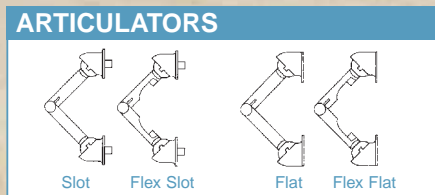


ARTICULATORS, TRAYS & SPINES

- Regular slot articulator..... PC 100010
- Flex slot articulator..... PC 100011
- Regular flat articulator..... PC 100013
- Flex flat articulator..... PC 100012
- Tray and Spine..... PC 100025

HINGES OR SOCKETS ONLY

- Regular Hinges Only.....PC 100020
- Flex Hinges Only.....PC 100021
- Flat Sockets Only.....PC 100023
- Slot Sockets Only.....PC 100022

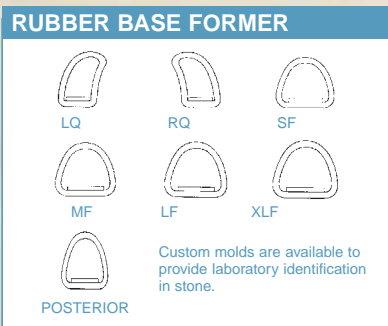


VERTEX RUBBER BASE FORMERS

- (Sold individually, mix or match)
- LQ, 16 MM (Regular).....PC 100050
 - LQ, 20 MM (Deep)..... PC 100060
 - RQ, 16 MM (Regular)..... PC 100051
 - RQ, 20 MM (Deep).....PC 100061
 - SF, 16 MM (Regular).....PC 100052
 - SF, 20 MM (Deep).....PC 100062
 - MF, 16 MM (Regular).....PC 100053
 - MF, 20 MM (Deep)..... PC 100063
 - LF, 16 MM (Regular).....PC 100054
 - LF, 20 MM (Deep).....PC 100064
 - XLF, 16 MM (Regular).....PC 100055
 - XLF, 20 MM (Deep).....PC 100065
 - Posterior.....PC 100056




VERTEX ACCESSORIES

- Articulator Adhesive, Super Thin 20 g.....PC 100030
- Articulator Adhesive, Extra Thick 1 oz.....PC 100033
- Accelerator, 2 fl. oz., with pump sprayer.....PC 100040
- Accelerator, 8 fl. oz., with pump sprayer.....PC 100041
- Pump Sprayer (Fits both Accelerator bottles).....PC 100090
- Round Carbide Bur HP.....PC 100070
- Plastic Stop Rod, Pkg of 20.....PC 100091
- Release Block.....PC 100044

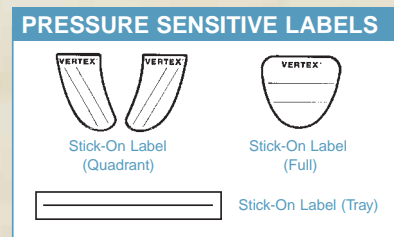

LABELS

- Stick-On Label (Left & Right Quadrant).....PC 100092
- Stick-On Label (Full Quadrant)..... PC 100093
- Stick-On Label (Tray).....PC 100094

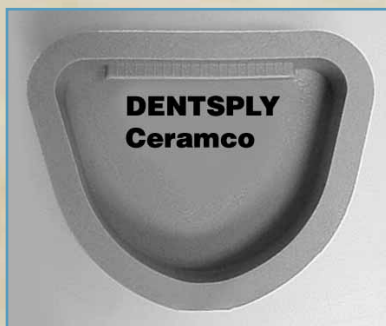

CUSTOM MOLDS

Custom molds are available to all labs. It provides your lab with the opportunity to promote your laboratory, while providing you with a professional and clean mold for your dentists.



DENTSPLY International
570 West College Ave., York, PA 17405
800-487-0100
www.ceramco.com or www.dentsply.com

Distirbuted by: Dentsply Canada,
161 Vinyl Ct., Woodbridge, Ontario L4L 4A3
905-851-5374