




Easy to fabricate, superior quality relines.

DuaLine® dual cure reline material makes quality in-office relines easy and patient pleasing. Its self-curing properties ensure an extremely accurate reline. And when cured in a DENTSPLY visible light curing unit, DuaLine reline material achieves its full physical properties and superior strength. Patients will appreciate the convenience of a single visit, on-the-spot reline and improved comfort of a better fitting denture.

Triad® DuaLine® dual-cure reline material makes it easy to achieve superior clinical quality relines

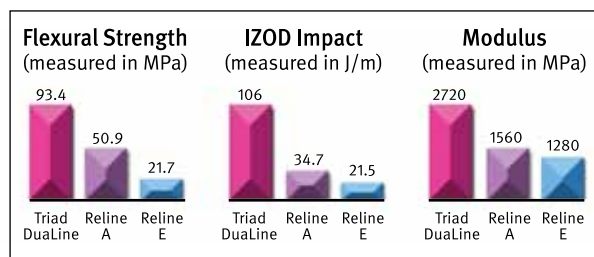

A stronger, better fitting reline

DuaLine dual-cure material is a breakthrough extension of the Triad material line, chemically engineered to be stronger than products containing MMA. DuaLine reline material is also stronger than other alternatives because it's less porous, and its full physical properties are efficiently achieved through post curing in a Triad® 2000 visible curing unit or Enterra® VLC curing unit. Its self-curing properties provide an extremely accurate impression, without distortion, ensuring a great fit.

No MMA or BMA



Like all Triad materials, DuaLine dual cure reline material is doctor and patient friendly. Because it contains no potentially sensitizing methyl methacrylate monomer or butyl methacrylate monomer, its virtually odor- free which is sure to please patients.

Eliminate porosity

DuaLine reline material isn't hand mixed. It's a paste/paste mixture that's extruded directly onto the denture. Since there's no air entrapment, there's virtually no staining, odor, or weakening of the reline.

Everybody's happy

DuaLine in-office relines are a business-building, patient-pleasing service. Reline patients' ill fitting, uncomfortable dentures right on-the-spot so they can leave with a better fitting denture after one short visit.



Part numbers:

- 905530 Intro Package
2 – 56g cartridges,
8 mixing tips,
20ml bonding agent,
15ml Air Barrier Coating,
25 brushes
- 905537 Refill package
2- 56g cartridges, 8 mixing tips
- 905539 Cartridge dispenser
1 dispensing gun, 1 plunger
- 70210 Triad® 2000 Visible Light Curing Unit
- 9495100 Enterra® VLC Unit
- 905538 Mixing Tip Refill
16 mixing tips
- 905220 Dispenser Plunger Refill
2 plungers
- 89699 Disposal Brushes
25 brushes

	Triad® 2000	Eclipse® EPU 2	Enterra® VLC Unit	interra™
Triad® Appliances				
Denture Baseplates/Record Bases	X	X	X	X
Impression Trays	X	X	X	X
Orthodontic Appliances	X	X	X	X
Denture Relines	X	X	X	X
Nightguard/Splint	X	X	X	X
Provisionals	X	X		
Full/Partial Denture	X	X		
Implant Placement Appliances	X	X	X	X
Name ID Tags	X	X		

DENTSPLY

Trubyte®
Triad® DuaLine®
Dual Cure Reline Material

DENTSPLY International, DENTSPLY Trubyte, 570 West College Ave., York PA. 17405-0872
Distributed By: DENTSPLY Canada, 161 Vinyl Court, Woodbridge, Ontario L4L 4A3
prosthetics.dentsply.com or dentsply.com