



FORM

Remarkable pen-style design that feels balanced in the hand and beautiful in its details.



Multifunctional charging base

The SmartLite® Pro comes with a modern-design charging base with a built-in Radiometer and room for extra tips. The intuitive battery management system includes two quick connect batteries for constant availability. The latest generation LiFePO4 battery technology ensures a full day of clinical operation with only one charge.

Ordering Information

644400	SmartLite® Pro Introductory Kit	1 x Handpiece 1 x Cure Tip 1 x Transillumination Tip 2 x Batteries 1 x Charging Base Accessories
644409	SmartLite® Pro Basic Kit	1 x Handpiece 1 x Cure Tip 2 x Batteries 1 x Charging Base Accessories
644401	SmartLite® Pro Battery Refill	
644402	SmartLite® Pro Sleeves Refill	
644403	SmartLite® Pro Eye Protection Shields Refill	
644404	SmartLite® Pro Power Connector Refill	
644405	SmartLite® Pro Transillumination Tip Refill	
644406	SmartLite® Pro PolyCure Tip Refill	
644407	SmartLite® Pro Cure Tip Refill	

NEW! SmartLite® Pro
Modular LED Curing Light

For more information visit www.dentsplysirona.com.
Find us on www.facebook.com/dentsplysirona.restorative

ML44002B (4/30/19)
© Dentsply Sirona 05/2019


THE DENTAL
SOLUTIONS
COMPANY™

**Dentsply
Sirona**



CLINICAL
PERFORMANCE
NEVER LOOKED
SO GOOD.

NEW! SmartLite® Pro
Modular LED Curing Light

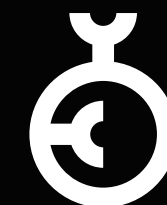

GERMAN
DESIGN
AWARD
NOMINEE
2020

**Dentsply
Sirona**

More than just a curing light

We weren't satisfied with creating one of the most advanced curing lights ever. It had to be one of the best-designed as well. The SmartLite® Pro curing light is a unique modular device in a remarkable, all-metal housing. But it delivers far more than good looks. State of the art optics provide leading quality of cure for reliable outcomes. SmartLite® Pro. **Designed to perform.**



⊕ Versatility

The SmartLite® Pro features a modular design with interchangeable tips for a variety of clinical indications. These quick connect, 360-degree rotational tips are easy to change. The forward-thinking design also allows for the addition of new tips as technology advances in the future.



⊕ Single-button simplicity

Intuitive one-button operation with audible and tactile signals allows you to focus on what really matters.

⊕ Metal housing

Medical-grade stainless steel and anodized aluminum provide robust durability, year after year.

⊕ Transillumination

The additional Transillumination Tip is a diagnostic aid for the visualization of interproximal caries, evaluation of cracked teeth and the illumination of endodontic access. It's sophisticated design includes a 1 mm focused light output with 2 different settings (anterior and posterior).


FUNCTION

Forward-thinking composition of clinical performance and modular versatility.

⊕ State of the art optics


Optimized 4 LED design provides a uniform light distribution over the whole active curing area. Irradiance: ~1.250 mW/cm²

SmartLite® Pro beam profile at 0 mm distance (Blue Light Range)¹

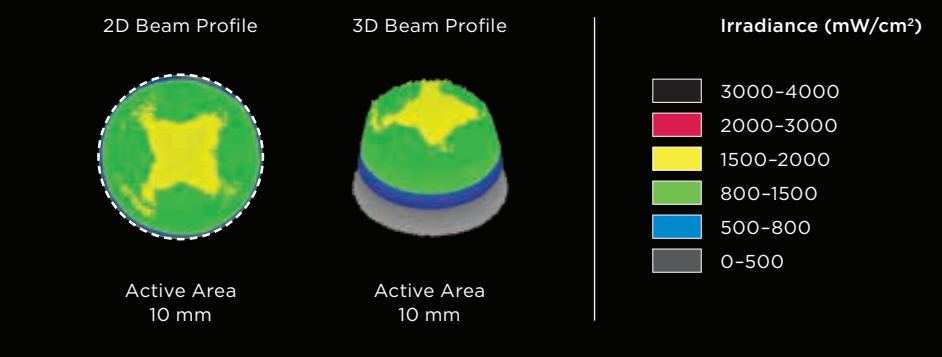

⊕ Reliable clinical outcomes


Large 10 mm active curing diameter to cover bigger restorations combined with an excellent beam collimation for reliable curing even over larger distances.

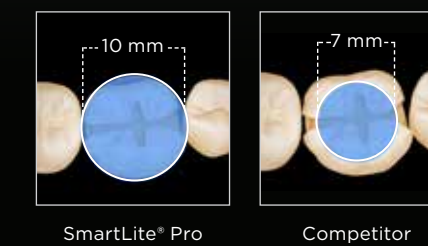

SmartLite® Pro: Collimated beam for reliable curing over larger distances.



Before

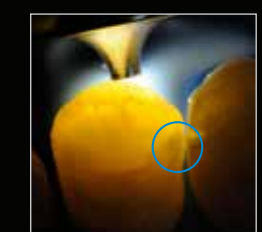

After