

+
**WE
KNOW
ENDO.**



WaveOne® Gold

Surf the canal
with confidence



Engineered with the latest reciprocating technology

The single-file simplicity of WaveOne® is now enhanced in WaveOne® Gold. Made with advanced metallurgy, it is the next generation of performance from Dentsply Sirona.



- Increases patient safety,* while reducing shaping time across a wide range of cases**
- Optimized with enhanced design features by leading endodontists
- Part of a complete solution that includes obturation products designed to match the shapes created by WaveOne® Gold files
- The WaveOne® Gold family now includes a single glide path file with reciprocating motion for simplicity from the start


Optimized tip

Variable taper

GOLD thermal treatment

Evolving cross-section

Reciprocating motion

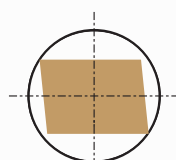

ISO color coded

variable taper



020 7% Size Small



025 7% Size Primary



035 6% Size Medium



045 5% Size Large

*compared to WaveOne®

**internal laboratory testing results

Simplicity for you, safety for patients

Reinforces patient safety

- Primary WaveOne® Gold file is 50% more resistant to cyclic fatigue than WaveOne® Primary file
- Reduced screwing effect** compared to standard rotary systems


Covers a wider* range of canal morphologies

- Enhanced file flexibility thanks to the Gold technology
- Extended size range (Small, Primary, Medium, Large)


Shortens the shaping time

- Takes the cutting efficiency to a higher level
- A single file per treatment also translates into time savings for canal shaping and irrigation


Cyclic fatigue resistance**

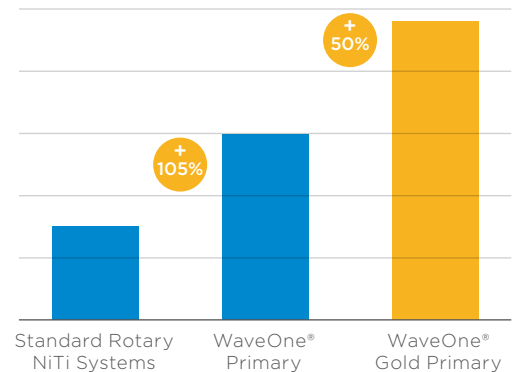

Flexibility**

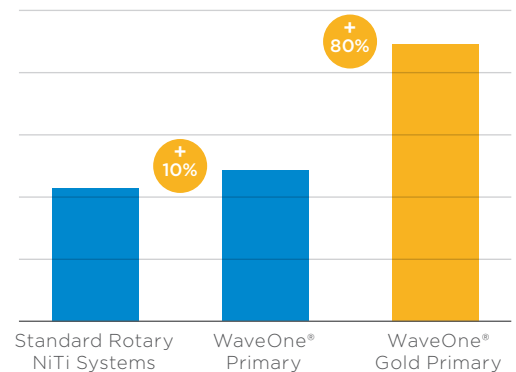

Shaping working time**

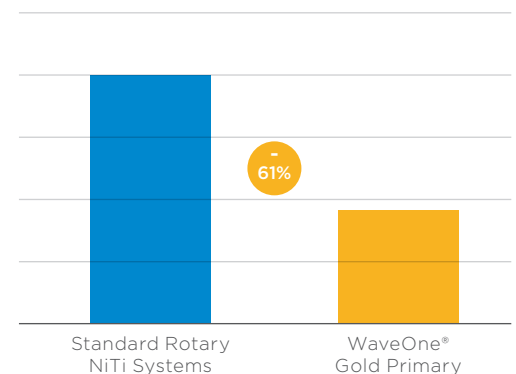

*compared to WaveOne®

**internal laboratory testing results



The next generation of performance

Dentsply Sirona has teamed up with the best endodontic professionals to develop WaveOne® Gold.



Courtesy of Dr. J. Webber

The design team speaks – Dr. C. Ruddle, Dr. J. Webber, Dr. W. Pertot, Dr. S. Kuttler

Shaping canals has never been easier

Since its launch in 2011, Dentsply Sirona's WaveOne® has set a new standard in root canal shaping. All around the world general dentists and endodontists have voiced their appreciation of the simplicity WaveOne® brings to their efforts to provide quality endodontic treatment to their patients.

We have created a new and exciting generation of WaveOne® to bring you WaveOne® Gold.

WaveOne® Gold brings to the endodontic procedure a metallurgically advanced single-file technique for shaping canals. Further, we optimized the tip diameters, tapers, and the cross-section to produce a file that really improves safety, efficiency and flexibility when preparing canals. WaveOne® Gold takes simplicity and treatment success to another level.

One complete solution

The WaveOne® Gold family of products provides a comprehensive treatment solution to promote confidence and predictable shape.

Motors

X-Smart® Plus and X-Smart IQ™ are pre-programmed for use with WaveOne® Gold for simplicity

Glide Path

The WaveOne® Gold Glider optimizes your glide path before shaping with WaveOne® Gold Primary**

Shaping

WaveOne® Gold shaping files are available in four sizes with three working lengths (21, 25, 31 mm) to treat a wide range of cases

Obturation

Dedicated WaveOne® Gold obturation products include paper points, variable tapered Conform Fit™ gutta-percha points and GuttaCore® Crosslinked Gutta-Percha Core Obturators

Dentsply Sirona Restorative Solutions

Please consult our website for restorative solutions
www.dentsplysirona.com

WaveOne® Gold Glider

Create your path with confidence




This new single glide path file has been designed in collaboration with a team of endodontist experts : Dr. E. Berrutti, Dr. G. Cantatore, Dr. A. Castellucci, Dr. S. Kuttler, Dr. W. Pertot, Dr. C. Ruddell, Dr. J. Webber, Dr. J. West.

Optimize your shaping preparation with a single reciprocating glide path file that is part of the WaveOne® Gold solution.


- Respects the natural shape of the canal
- Increases patient safety with enhanced flexibility and resistance to cyclic fatigue
- Sterile and ready for single use, which preserves cutting efficiency, reduces breakage and prevents cross contamination

The WaveOne® Gold Glider is only for use in reciprocating motion with compatible Dentsply Sirona motors. Helpful pre-programmed motor settings are available.

Lower bending resistance than PathFile, for an increased flexibility thanks to the gold heat treatment**

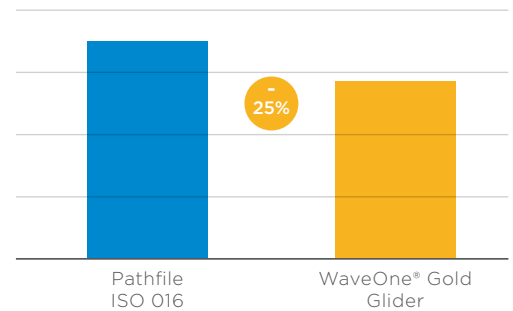

Better resistance to cycle fatigue**

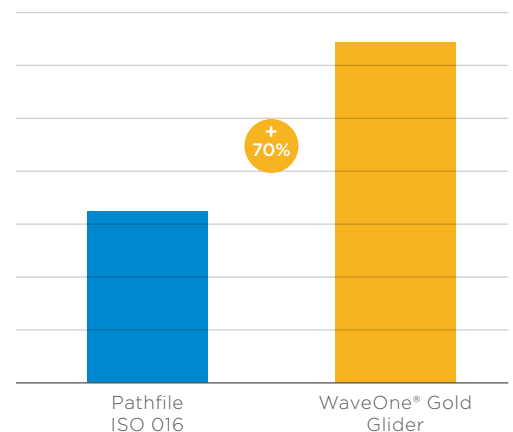

One sequence for most of your cases

Scouting

ReadySteel®



010 02 - 25mm

Glide Path

WaveOne® Gold Glider



015 02 - 25mm
Variable Taper


Shaping

WaveOne® Gold



Size Primary,
025 07 - 25mm
Variable Taper

Same cross section as the WaveOne® Gold shaping files

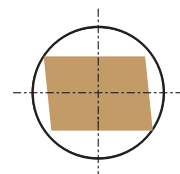

**internal laboratory testing results

Expanded selection of corresponding products

With the addition of a glide path file, WaveOne® Gold gives you the convenience of an entire family of products designed to work together, right from the start.

Original WaveOne® Gold family of products



Expanded WaveOne® Gold family of products under our new Dentsply Sirona branding



The new gold standard of gutta-percha introduced for WaveOne® Gold




Conform Fit™

WaveOne® Gold Gutta-Percha Master Cones are latex-free and reformulated to match the precise shapes created by the WaveOne® Gold instruments.

- Improved heat transfer up to 4mm
- Variable taper
- Precise fit and accurate tug back
- Latex free for patient safety
- Part of a system based solution

Superior apical gutta-percha point filling volume**

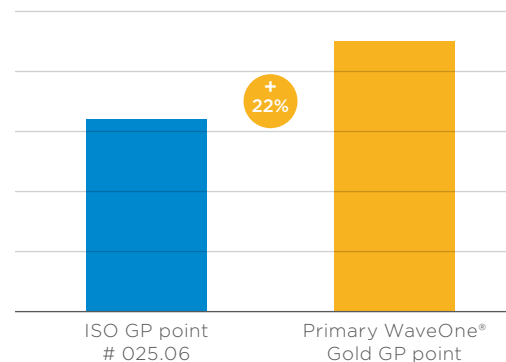

Clearly identifiable: same color code as the file. to traditional hand rolled gutta-percha with the



Paper points and GuttaCore® for WaveOne® Gold



Paper points

- Paper points are packaged in individual sterile packs of 5 points
- Ready to use, part of a total system solution for simple identification




GuttaCore® for WaveOne® Gold

GuttaCore® obturators feature a crosslinked gutta-percha core that sends the surrounding warm, flowable gutta-percha 3-dimensionally to follow the natural canal anatomy.

- An efficient technique that offers the ease of a single insertion, so your patients spend less time in your chair.
- The flowability of warm gutta-percha to achieve a complete fill throughout the canal system.
- The confidence of knowing GuttaCore® travels where traditional cold gutta-percha stops short.

4 sizes (small, primary, medium, large) and improved dimensional stability up to 40% compared related WaveOne® Gold files.

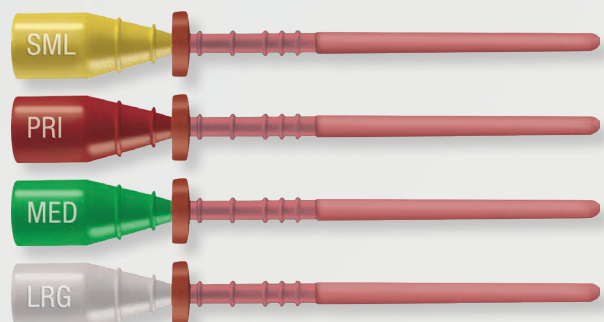



Product references

Reciprocating glide path files



Glider - 015, 02 variable taper

3x	21mm	B ST WIGG 321 015
3x	25mm	B ST WIGG 325 015
3x	31mm	B ST WIGG 331 015
6x	21mm	B ST WIGG 621 015
6x	25mm	B ST WIGG 625 015
6x	31mm	B ST WIGG 631 015

Reciprocation shaping files



Small - 020, 07 variable taper

3x	21mm	A 0753 221 GOS 03
3x	25mm	A 0753 225 GOS 03
3x	31mm	A 0753 231 GOS 03
6x	21mm	A 0756 221 GOS 03
6x	25mm	A 0756 225 GOS 03
6x	31mm	A 0756 231 GOS 03

Primary - 025, 07 variable taper

3x	21mm	A 0753 221 GOP 03
3x	25mm	A 0753 225 GOP 03
3x	31mm	A 0753 231 GOP 03
6x	21mm	A 0756 221 GOP 03
6x	25mm	A 0756 225 GOP 03
6x	31mm	A 0756 231 GOP 03

Medium - 035, 06 variable taper

3x	21mm	A 0753 221 GOM 03
3x	25mm	A 0753 225 GOM 03
3x	31mm	A 0753 231 GOM 03
6x	21mm	A 0756 221 GOM 03
6x	25mm	A 0756 225 GOM 03
6x	31mm	A 0756 231 GOM 03

Large - 045, 05 variable taper

3x	21mm	A 0753 221 GOL 03
3x	25mm	A 0753 225 GOL 03
3x	31mm	A 0753 231 GOL 03
6x	21mm	A 0756 221 GOL 03
6x	25mm	A 0756 225 GOL 03
6x	31mm	A 0756 231 GOL 03

Assortment Reciprocating File 1x Small/Primary/Medium/Large

21mm	A 0754 221 G9A 03
25mm	A 0754 225 G9A 03
31mm	A 0754 231 G9A 03

WaveOne® Gold blister sequence

25mm	B ST WIGS 325 AST
------	-------------------

WaveOne® Gold absorbent points



Box of 180x Sterile Paper Points

Small	A 175W 000 00S 03
Primary	A 175W 000 00P 03
Medium	A 175W 000 00M 03
Large	A 175W 000 00L 03

Conform Fit™ gutta-percha points for WaveOne® Gold



Box of 60x Variable Tapered Gutta Percha Points

Small	A 175X 000 00S 00
Primary	A 175X 000 00P 00
Medium	A 175X 000 00M 00
Large	A 175X 000 00L 00

GuttaCore® for WaveOne® Gold



Blister of 6x Obturators

Small	B NS GCWG 500 SML
Primary	B NS GCWG 500 PRI
Medium	B NS GCWG 500 MED
Large	B NS GCWG 500 LRG

Box of 30x Obturators (5 Blisters)

Small	B NS GCWG C00 SML
Primary	B NS GCWG C00 PRI
Medium	B NS GCWG C00 MED
Large	B NS GCWG C00 LRG

Kits and motors

X-SMART IQ WaveOne Gold Starter Kit (W1GG)

B 00 XSIQ WIG KIT

- 1 Motor with Contra-Angle
- 1 Pack of 15 Primary WaveOne Gold Files
- 1 Pack of 16 Assorted WaveOne Gold Files
- 1 Blister of 3 WaveOne Gold Glider
- 1 iPad Protective Case

Motors

X-SMART IQ™ Motor	A1054
X-SMART® Plus Motor	A1032

Oven ThermaPrep® 2



230V (EU)	A 0176 230 002 00
230V (UK)	A 0176 230 002 33
115V	A 0176 115 002 00