

# Atlantis® CustomBase Solution— Together we make implants personal



“Ordering the abutment and crown, designed together, in a multilayer material, gives us the foundation we need to finalize the restoration just like natural teeth.”

Kungskronan Dental, Gothenburg, Sweden

## Now available with cut-back options and Cercon® xt ML

The premium quality Atlantis CustomBase Solution combines a patient-specific abutment and crown for screw-retained restorations. Simultaneous design and manufacturing of the abutment and crown ensure best possible results. This gives you peace of mind and allows you to focus on what you do best—finalizing the implant restoration.

**Cut-back options**—The Atlantis Crown, cut-back can be ordered with 4 different types of cut-back designs enabling case specific flexibility between esthetics and function.



**Cercon xt ML**—With the Atlantis Crown in Cercon xt ML multilayered material, esthetics are enhanced even further. Select the ideal gradient position to precisely match your patient's tooth shade and design natural esthetic results.



# Mimic nature without compromising function

Milled in multilayered zirconia disks, Cercon xt ML has extra-high translucency (approximately 49%) to suit patient dentition.

Atlantis Crown is strong and suitable for single-unit restorations in the anterior and posterior region (flexural strength: 750 MPa).

16 VITA shades plus 1 bleached shade with its own incisal shade for natural esthetic results.



## Cut-back options

- **Full cut-back**—For full porcelain veneering
- **Lingual collar**—For easier handling with pliers
- **Proximal contacts**—To minimize risks for porcelain cracks in the proximal area
- **Facial cut-back**—Allows for full contour zirconia in the occlusal and proximal areas and porcelain veneering in the esthetic area




Full cut-back

Lingual collar

Proximal contacts

Facial cut-back


Make your selection in Atlantis WebOrder!

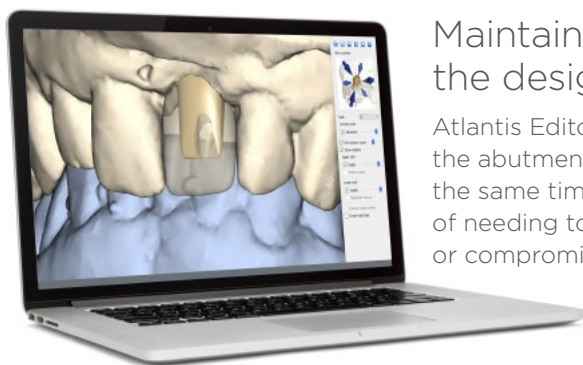

Atlantis Editor

## Maintain control throughout the design process

Atlantis Editor enables you to modify the abutment and crown design at the same time. This minimizes the risk of needing to alter the abutment design or compromise on the crown design.

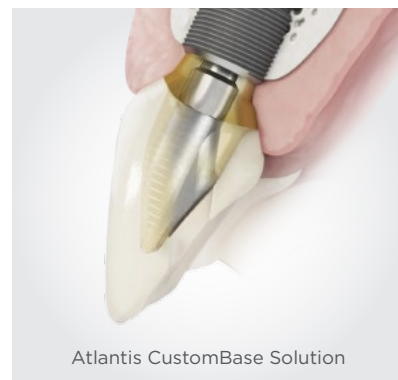

Atlantis CustomBase Solution



## Offer your patients superior esthetics

Please contact your sales representative or customer service team for more information about Atlantis CustomBase Solution and Atlantis Crown in Cercon xt ML material.

[www.dentsplysirona.com](http://www.dentsplysirona.com)

THE DENTAL  
SOLUTIONS  
COMPANY™

 **Dentsply  
Sirona**