Cleanup

Most cements give you limited time to clean up, increasing the likelihood of excess cement being left on the tooth. The Calibra® Cements offers easy cleanup to ensure removal of excess cement.

Selection

Determining which cement is best suited for a specific indication has been complicated and confusing, and can cause a restoration to fail if selected incorrectly. Cement selection doesn't have to be so hard.

Why spend so much time and effort on a restoration, only to have its longevity compromised by incomplete cement cleanup or incorrect selection?

Simple selection means successful restorations.

That's why we've made it easy to choose the right Calibra® cement for any indication. Relying on Calibra® cements, the trusted brand with over a decade of clinical performance, for all your cementation needs can also simplify ordering and inventory, as well as reduce the number of techniques and delivery systems your team must master. With the Calibra® family of definitive cements, selection is simple and success is inevitable.



Calibra® Universal Self-Adhesive Resin Cement

The routine, self-adhesive cement for nearly any indication²

Calibra® Ceram Adhesive Resin Cement

The maximum strength adhesive cement with immediate and high bonding for ceramics and CAD/ CAM restorations³

Calibra® Veneer Esthetic Resin Cement

The enduring, esthetic veneer cement with 12 years of proven clinical performance⁴

Calibra[®] Bio Bioceramic Luting Cement

Bioactive cement.

Self-repairing hydroxyapatite (HA) layer. For long lasting marginal integrity, naturally.

titter, Robert. "A Practical, Cost-effective Approach to Adhesive Cementation." Dental Economics. N.p., n.d. Web. metal-based inlays, onlays crowns, bridges, and posts.

interal-based intays, onlays crowns, bridges, and posts.

3. Internal data on file. For more information, contact Consumables-Data-Requests@dentsplysirona.com.

4. Farah, John W. and Ireland, Jennifer. "Twelve-year Performance of Calibra Esthetic Resin Cement". The Dental Advisior, number 27: May 2012.

Calibra® Universal Self-Adhesive Resin Cement

Dual Cure AutoMix Syringe Refill Package

2 Syringes (4.5g each), 20 Mixing Tips

607402 O Light

607403 Medium 607405 () Translucent 607406 Opaque

607407 () Bleach



607086 **Dual Cure AutoMix Syringe Mixing Tip Refill** (50)

Calibra® Ceram Adhesive Resin Cement

607100 **Combo Kit**

1 AutoMix Syringe (4.5g) - Translucent Shade 10 Mixing Tips

1 Bottle Prime&Bond elect® Adhesive (5mL)

25 Flocked Applicator Tips 1 Dispensing Well

Dual Cure AutoMix Syringe Refill Package

1 Syringe (4.5g), 10 Mixing Tips

607191 O Light

607192 Medium

607194 Translucent 607195 Opaque

607196 O Bleach

607086 **Dual Cure AutoMix Syringe Mixing Tip Refill** (50)

Calibra® Bio Bioceramic Luting Cement

607601 Introductory Kit

7 Capsules (334mg powder, 196mg liquid ea) 1 Capsule Extruder 2

THE DENTAL

SOLUTIONS

COMPANY

607604 Capsule Refill

14 Capsules (334mg powder, 196mg liquid ea)

606.06.596 **Capsule Extruder 2**

5546068 Dentsply Sirona Capsule Mixer

27700010000 PROMIX one™ Capsule Mixer

©2023 Dentsply Sirona Inc. All rights reserved ML070007B (R 1/27/23)

For more information and to request a demo, contact your Dentsply Sirona Representative.

Dentsply Sirona | 13320-B Ballantyne Corporate Place | Charlotte, NC 28277 | 1-800-532-2855 | www.dentsplysirona.com

Esthetic Resin Cement

Calibra® Veneer

5 Syringes Light Cure (2g each)

1 each Shade: Light, Medium, Translucent, Opaque, Bleach

5 Syringes Try-In Paste (1.8g each) 1 each Shade: Light, Medium, Translucent, Opaque, Bleach

1 Prime&Bond XP Adhesive (5mL)

1 Syringe Calibra[®] Silane Coupling Agent (3mL) 1 Syringe Caulk® Tooth Conditioner Gel (3mL)

25 Applicator Needles

50 Flocked Applicator Tips 1 Dispensing Well

1 Light Protection Cover

Refill Package

1 Syringe (2g)

607201 Light Shade

607202 Medium Shade 607204 Translucent Shade 607205 Opaque Shade

607206 Bleach Shade

607300 Try-In Paste Accessory Pack

5 Syringes Try-In Paste (1.8g each)

1 each Shade: Light, Medium, Translucent, Opaque, Bleach

Try-In Paste Refill Package

607306 O Bleach

607301 () Light 607302 Medium 607304 O Translucent 607305 Opaque



Calibra® Cements

The Simple Choice for Easy Cleanup



THE DENTAL

SOLUTIONS

COMPANY



Easy Cleanup and Simple Selection for Successful Restorations

Easy cleanup means successful restorations.

97% OF DENTISTS SAY CLEANING UP EXCESS CEMENT IS IMPORTANT TO ENSURE HEALTH OF THE GINGIVA1

Failure to clean excess cement can be one of the contributing factors² to crown failures and secondary caries. Calibra® cements offer easy, unhurried cleanup for convenience and the assurance of a long-lasting restoration.









With the easy cleanup of the Calibra cements family, you can make sure you're no longer one





Tack Cure with Confidence

A wide tack cure window of up to 10 seconds allows for the speed of tack curing without the worry of over-curing.^{3, 4}



A 45-second extended gel phase gives you the time you need for thorough cleanup of excess cement.4


We've Simplified Cement Selection

Use our Selection Tool to find the right cement and ensure the success of your restorations. For our comprehensive cementation selection guide, visit us at

https://www.dentsplysirona.com/en-us/discover/discover-by-brand/calibra-cements.html

Material	Indication	Preparation	Calibra® Cement
Composite/Glass Ceramic	Crown	Retentive	Calibra Ceram
		Non-Retentive	Calibra Ceram
High Strength Glass Ceramic	Crown, Bridge	Retentive	Calibra Universal
		Non-Retentive	Calibra Ceram
Zirconia/Metal	Crown, Bridge	Retentive	Calibra Universal
		Non-Retentive	Calibra Ceram

Endodontic Post	Calibra Universal
Inlay, Onlay	Calibra Ceram
Veneer	Calibra Veneer

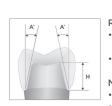

• Preparation with adequate cervical-occlusal height (>4mm)

• Taper between 10°-20°

Non-Retentive²

- Short clinical crown (≤4mm)
- Over-tapered preparation greater than 20°

Calibra® Dentsp Sirona

Universal

CALIBRA® UNIVERSAL SELF-ADHESIVE RESIN CEMENT THE EASY, UNIVERSAL CHOICE

With a formulation that suits so many types of restorations, it's no wonder it's called Calibra Universal cement. It's the simple choice for successful results across a wide array of indications.

- One-step cementation means no additional etchant or adhesive
- High immediate 6-minute bond strength for a strong initial bond²
- 5 esthetic shade options with Shade Stable® technology virtually eliminates color shift over time²

2. Internal data on file, For more information, contact Consumables-Data-Requests@dentsplysirona.com

CALIBRA® CERAM ADHESIVE RESIN CEMENT

THE MAXIMUM STRENGTH CHOICE

Non-retentive preparations, requires high retentive bond strength.

Calibra® Ceram cement offers strong bonds and long-term success in all-ceramics and CAD/CAM restorations.

- Simple system for superior clinical performance. No need for a self cure activator, and no mixing of primers when used with Prime&Bond elect®
- Immediate and long-term high retentive bond strength on ceramics compared to leading brands means successful restorations and greater
- Low film thickness

Shade Stability


• 5 esthetic shade options with Shade Stable® technology virtually eliminate color shift over time

Calibra Multilink® Variolink® RelvX™

Test: Color stability. All cements were light cured. Samples were stored in 37°C distilled water for 24hr (baseline) and 60°C distilled water for 8 weeks. Source: Dentsply Sirona internal data.²




Test: Shear Bond Strength. Calibra Ceram Cement was used with Prime&Bond elect Adhesive; other materials with respective prime and adhesives according to manufactures' instructions. All cements

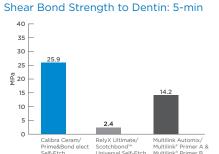

were self cured. Specimens were debonded after 5 minutes storage

Crown Retention After 10000 TC

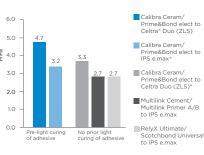

Test: Crown-pull Bond Strength Testing, Ceramics were treated wit 5% HF and silanated. All adhesives were used in self-etch mode. All cements were self-cured except for the group with optional light-curing through crown. Artificial aging: 10000 thermocycling

* light curing of adhesive through crown



CALIBRA® VENEER ESTHETIC RESIN CEMENT THE CLINICALLY3 PROVEN CHOICE

With its unbeatable track record, Calibra Veneer Cement is the only cement you'll ever need for successful, long-lasting veneer restorations. • Over a decade of proven clinical performance

- A visible light cured cement with virtually unlimited cleanup time²
- Thixotropic nature allows for favorable handling
- Try-in pastes closely match their corresponding shades
- 5 esthetic shade options with Shade Stable® technology virtually eliminate color shift over time

3. Farah, John W. and Ireland, Jennifer. "Twelve-year Performance of Calibra Esthetic Resin Cement". The Dental Advisior, number 27: May 2012.

 $The Key Group, Inc.\ 2015\ Dental Omnibus: Quarter\ 1\ Voice of Customer Reports. For more information, contact Consumables-Data-Requests @dentsplysirona.com and the Consumable of Customer Reports of Cust$

RelyX Ultimate, Multillink Automix, Variolink Esthetic, IPS e.max, Scotchbond Universal, and Multillink Primer A/B are not registered trademarks of Dentsply Sirona.

^{2.} Amani, L. C., Pereira, C. A., Guerra, E., Souza, R., Jorge, A., & Bottino, M. (2012). Morphology and bacterial colonisation of tooth/ceramic restoration interface after different cement excess removal techniques. Journal of Dentistry, 40, 742-749.

^{4.} For excess cement cleanup, monowave output LED lights with a single peak output around 470nm are recommended. High power, dual or broad spectrum lights may cause premature hardening of excess cement

Check curing light effect on mixed cement prior to clinical use