

Midwest Stylus® Plus

High-Speed Handpiece

THE NEXT GENERATION AIR TURBINE FROM DENTSPLY SIRONA

More Power

All new rotor design provides up to 30% more power than previous Stylus generations.

More Control

Dynamic speed control in the head helps maintain cutting power while also protecting against excessive speeds that can cause premature bearing wear.

Less Noise

Unique bearing design, computer balanced rotor, and tuned air flow work together to minimize vibration and noise for a more pleasant work environment.

Improved Design


Ceramic bearings are protected with quick-stop technology to help keep a clean service record; 2-year warranty.



Fiber Optic | Push Button Chuck | Multi-port Water Spray


- 792400 Midwest Stylus® Plus — Midwest backend
- 792410 Midwest Stylus® Plus — Sirona backend
- 792420 Midwest Stylus® Plus — NSK backend
- 792430 Midwest Stylus® Plus — K-style backend

Call your Midwest Specialist at 800-800-2888 to schedule a demo today. Ask about special trade-in offers!

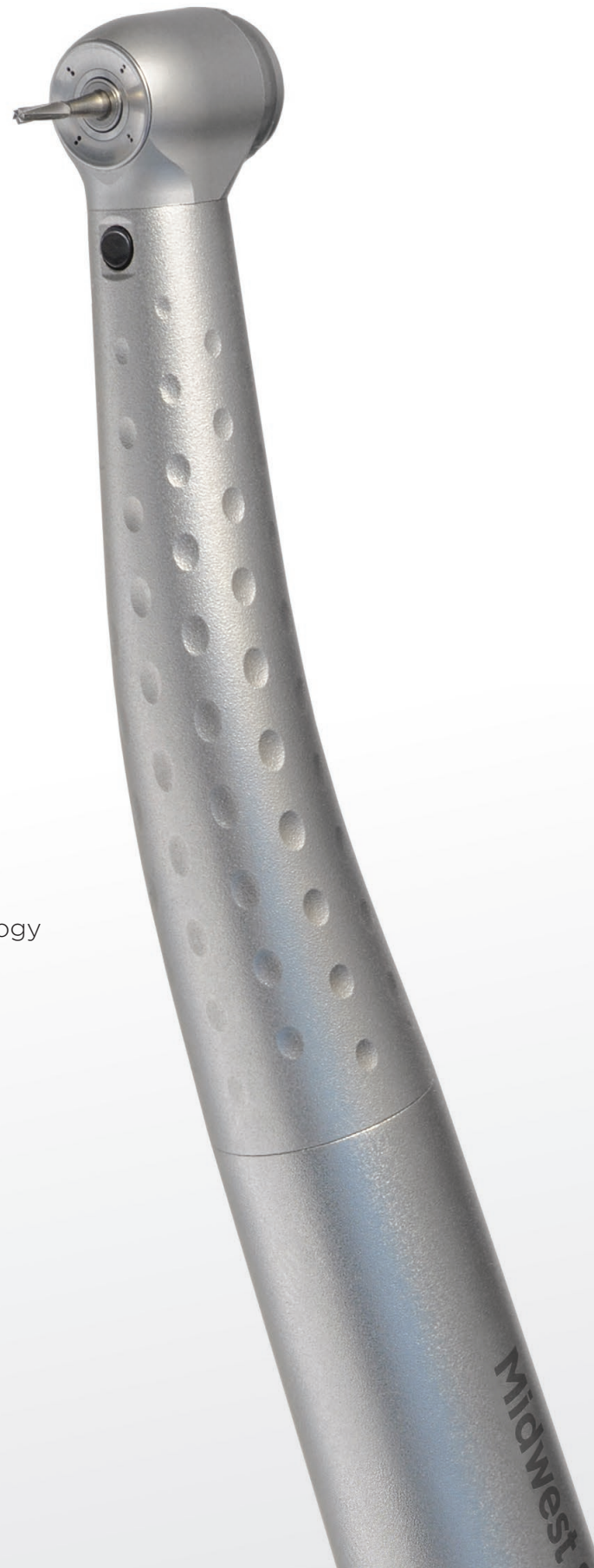

PERFORMANCE UNDER PRESSURE

Stay in control with the Midwest Stylus® Plus. Dynamic Speed Control in the head limits free speed to protect bearings against excessive RPMs at higher pressures. This also translates to lower noise without impacting peak cutting power.

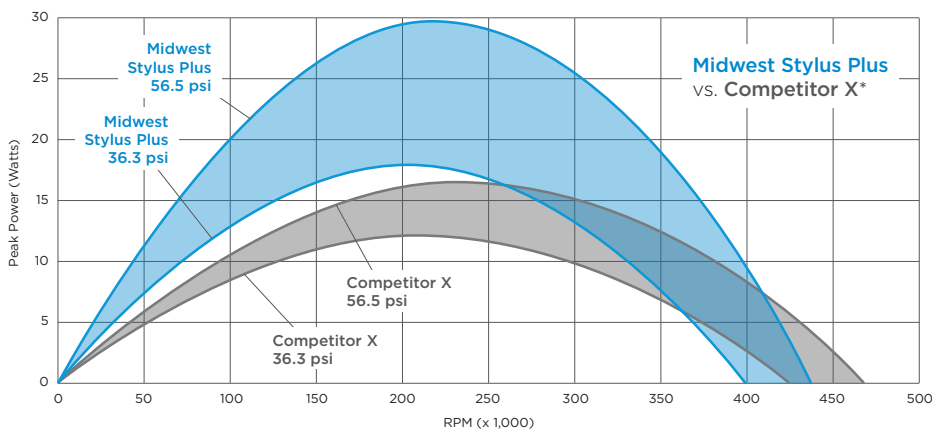

COMPLETE CUTTING SOLUTIONS

As a leader in both handpiece and bur manufacturing, Midwest offers the ultimate integration of design, performance, and integrity. Explore our full array of carbide and diamond burs with the [Midwest Bur Advisor](#), a virtual bur wallet and utilities program available for Apple, Android and Windows. Download it free at the Apple Store, on Google Play, or through www.MidwestBurApp.com.



KEEP IT CLEAN

In clinical studies, correct cleaning and lubrication have proven to be the most important factors for handpiece performance and longevity. To ensure you get the most out of your handpieces, the Midwest® Automate® system:

- Reduces handling time by more than 80%
- Automatically removes debris, precisely lubricates all moving parts, and accurately expels excess
- Ensures consistent maintenance every time
- Ends the excess mess of manual handpiece maintenance
- Optimizes long-term performance with quality maintenance



ORDER INFORMATION

Midwest Stylus® Plus High-Speed Handpieces

- 792400 Midwest Stylus® Plus — Midwest backend
- 792410 Midwest Stylus® Plus — Sirona backend
- 792420 Midwest Stylus® Plus — NSK backend
- 792430 Midwest Stylus® Plus — K-style backend

Couplers

- 790245 5-Hole Coupler — Midwest style
- 790246 6-Pin Coupler — Midwest style
- 790335 KS5 5-Hole Coupler — K-style
- 790336 KS6 6-Pin LED Coupler — K-style

Midwest® Automate® Automated Maintenance Unit

- 380200 Midwest® Automate® Unit
- 380090 Midwest Stylus® Maintenance Adapter
- 380450 K-style Maintenance Adapter
- 380480 NSK Maintenance Adapter
- 380510 Sirona Maintenance Adapter
- 380520 Absorption Pad (6 Pack)
- 380080 Midwest Plus™ Aerosol Spray Can

Warranty

Midwest Stylus® Plus is backed by nearly 90 years of Midwest handpiece experience, and a 2-year warranty.