

TPH Spectra® ST Effects

Universal Composite Restorative


Next-Level Esthetics, Simplified

For esthetically demanding anterior cases, TPH Spectra® ST Effects composite works with TPH Spectra® ST composite in a simplified layering technique.

- TPH Spectra® ST Effects opaque **dentin** shades are paired with TPH Spectra® ST **universal CLOUD shades** when a case calls for added opacity.
- One translucent **enamel** shade creates a highly translucent incisal edge when needed.
- Simplified layering. Made possible by **SphereTEC® technology**.


1 out of 3 people say that the first part of someone's face they notice is their teeth.²

OVER 1/3 of clinicians consider color matching or overall esthetics to be the most challenging aspect of an anterior composite restoration.³

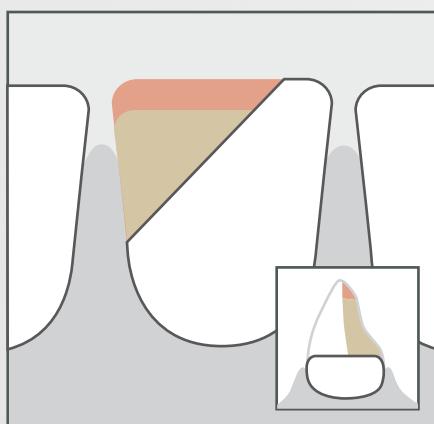


Opaque Dentin Shade

Restore lingual 1/2 of tooth thickness with opaque dentin and universal CLOUD shade.

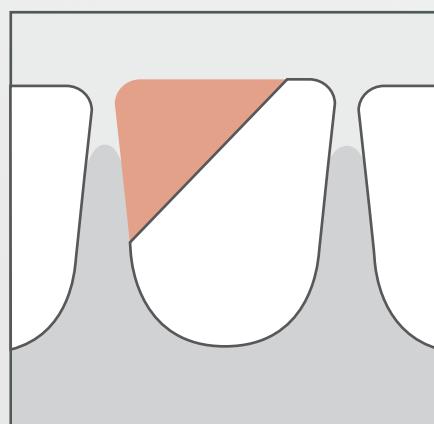

Universal CLOUD Shade

Restore facial 1/2 of tooth thickness with universal CLOUD shade.

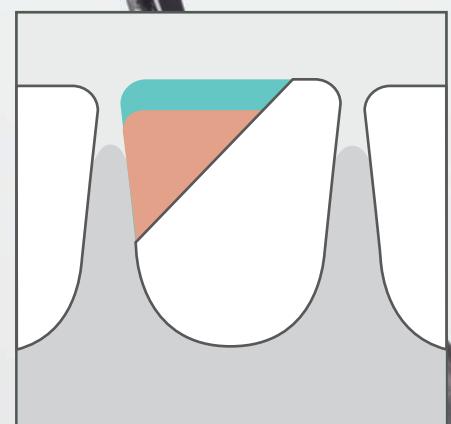

Translucent Enamel Shade (Optional)

For a highly translucent incisal edge, optionally use the translucent enamel shade for the facial incisal edge increment.

1. VITA is not a registered trademark of Dentsply Sirona.


2. <http://www.melbourneorthodontics.com/educational/study-7-ways-smiles-affect-first-impressions/>

3. 2018 The Key Group, Inc.

TPH Spectra® ST Effects

Universal Composite Restorative

Streamline your composite inventory with
TPH Spectra® ST + TPH Spectra® ST Effects Simplified Layering Technique.

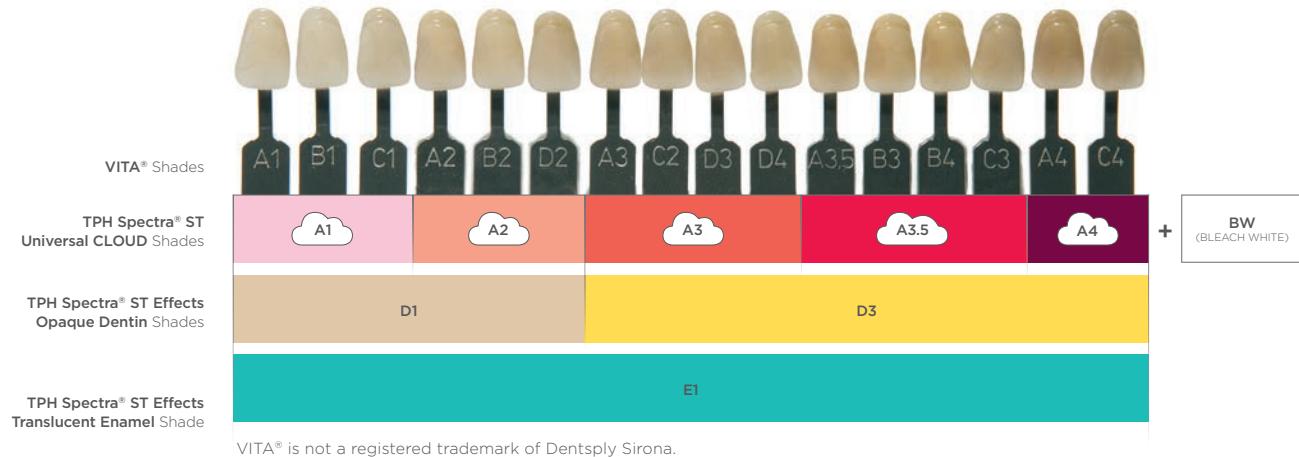

Ordering Information

TPH Spectra® ST Effects Universal Composite Restorative - Compules® Tips - HIGH VISCOSITY

642542 **Compules® Tips Introductory Kit**
28 Compules® Tips (0.25g each)
12 Shade D1
8 Shades D3
8 Shade E1
Illustrated Technique Guide
i-Shade Label



TPH Spectra® ST Effects Universal Composite Restorative - HIGH VISCOSITY REFILLS

HV Refills	Compules® Tips	Syringe
Shades	(16 x 0.25g)	(1 x 3g)
D1 (use with A1; A2)	642543	642546
D3 (use with A3; A3.5; A4)	642544	642547
E1 (use with A1; A2; A3; A3.5; A4)	642545	642548



642549 **TPH Spectra® ST Shade Guide Expansion Pack**
3 Shade Tabs (D1, D3, E1)



For more information and to learn the technique, contact your Dentsply Sirona Restorative specialist, or visit www.dentsplysirona.com/TPHSpectraSTEffects

¹ VITA is not a registered trademark of Dentsply Sirona.